

(12) INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(19) World Intellectual Property  
Organization  
International Bureau



(43) International Publication Date  
10 September 2004 (10.09.2004)

PCT

(10) International Publication Number  
**WO 2004/076485 A1**

(51) International Patent Classification<sup>7</sup>: **C07K 16/06**,  
1/36 // 17/02

PH, PL, PT, RO, RU, SC, SD, SE, SG, SK, SL, SY, TJ, TM,  
TN, TR, TT, TZ, UA, UG, US, UZ, VC, VN, YU, ZA, ZM,  
ZW.

(21) International Application Number:  
PCT/EP2004/002041

(84) Designated States (unless otherwise indicated, for every  
kind of regional protection available): ARIPO (BW, GH,  
GM, KE, LS, MW, MZ, SD, SL, SZ, TZ, UG, ZM, ZW),  
Eurasian (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM), Euro-  
pean (AT, BE, BG, CH, CY, CZ, DE, DK, EE, ES, FI, FR,  
GB, GR, HU, IE, IT, LU, MC, NL, PL, PT, RO, SE, SI, SK,  
TR), OAPI (BF, BJ, CF, CG, CI, CM, GA, GN, GQ, GW,  
ML, MR, NE, SN, TD, TG).

(22) International Filing Date: 1 March 2004 (01.03.2004)

(25) Filing Language: English

(26) Publication Language: English

(30) Priority Data:  
0304576.2 28 February 2003 (28.02.2003) GB  
60/464973 24 April 2003 (24.04.2003) US

**Declaration under Rule 4.17:**

— as to applicant's entitlement to apply for and be granted  
a patent (Rule 4.17(ii)) for the following designations AE,  
AG, AL, AM, AT, AU, AZ, BA, BB, BG, BR, BW, BY, BZ,  
CA, CH, CN, CO, CR, CU, CZ, DE, DK, DM, DZ, EC, EE,  
EG, ES, FI, GB, GD, GE, GH, GM, HR, HU, ID, IL, IN, IS,  
JP, KE, KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MA,  
MD, MG, MK, MN, MW, MX, MZ, NA, NI, NO, NZ, OM,  
PG, PH, PL, PT, RO, RU, SC, SD, SE, SG, SK, SL, SY, TJ,  
TM, TN, TR, TT, TZ, UA, UG, UZ, VC, VN, YU, ZA, ZM,  
ZW, ARIPO patent (BW, GH, GM, KE, LS, MW, MZ, SD,  
SL, SZ, TZ, UG, ZM, ZW), Eurasian patent (AM, AZ, BY,  
KG, KZ, MD, RU, TJ, TM), European patent (AT, BE, BG,  
CH, CY, CZ, DE, DK, EE, ES, FI, FR, GB, GR, HU, IE, IT,  
LU, MC, NL, PL, PT, RO, SE, SI, SK, TR), OAPI patent (BF,  
BJ, CF, CG, CI, CM, GA, GN, GQ, GW, ML, MR, NE, SN,  
TD, TG)

(71) Applicant (for all designated States except US): **LONZA  
BIOLOGICS PLC**, [GB/GB]; 228 Bath Road, Slough SL  
14DY, Berkshire (GB).

(72) Inventors; and

(75) Inventors/Applicants (for US only): **BONNERJEA**,  
Julian [GB/GB]; 4 South Maundin, Hughenden Valley,  
High Wycombe, Bucks HP14 4LZ (GB). **PRENETA**,  
Anna [GB/GB]; 35 Walton Park, Walton-on-Thames,  
Surrey KT12 3EU (GB).

(74) Common Representative: **LONZA BIOLOGICS PLC**;  
Legal Department, 228 Bath Road, Slough SL 14DY, Berk-  
shire (GB).

**Published:**

— with international search report  
— before the expiration of the time limit for amending the  
claims and to be republished in the event of receipt of  
amendments

For two-letter codes and other abbreviations, refer to the "Guid-  
ance Notes on Codes and Abbreviations" appearing at the begin-  
ning of each regular issue of the PCT Gazette.

(81) Designated States (unless otherwise indicated, for every  
kind of national protection available): AE, AG, AL, AM,  
AT, AU, AZ, BA, BB, BG, BR, BW, BY, BZ, CA, CH, CN,  
CO, CR, CU, CZ, DE, DK, DM, DZ, EC, EE, EG, ES, FI,  
GB, GD, GE, GH, GM, HR, HU, ID, IL, IN, IS, JP, KE,  
KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MA, MD,  
MG, MK, MN, MW, MX, MZ, NA, NI, NO, NZ, OM, PG,

(54) Title: ANTIBODY PURIFICATION BY PROTEIN A AND ION EXCHANGE CHROMATOGRAPHY

(57) Abstract: A novel method for selectively removing leaked protein A from antibody purified by means of protein A affinity chromatography is disclosed.



WO 2004/076485 A1

## ANTIBODY PURIFICATION BY PROTEIN A AND ION EXCHANGE CHROMATOGRAPHY

The present invention relates to the field of antibody purification in biotechnological production. It is an object of the present invention to describe a novel process for  
5 purification of such antibody.

Protein A chromatography is widely used in industrial manufacturing of antibodies since allowing for almost complete purification of antibodies, that is usually IgG, in a single step from cell culture supernatants. Protein A affinity columns inevitably are subject to some  
10 degree of leakage of ligand from the column upon repeated runs. Partly, this may be due to proteolytic clipping of protein A from the column; in industrial manufacture of antibody for pharmaceutical applications, no protease inhibitor cocktails may be added for regulatory reasons. Unfortunately, this protein A or protein A fragment contaminants retain their affinity for IgG and are difficult to remove from the purified antibody due to ongoing  
15 complex formation. Removal is mandatory since protein A which is a bacterial protein will elicit an unwanted immune response; further, model complexes formed by adding protein A to monomeric IgG have been reported to activate Fc-bearing leukocytes and the complement system to generate oxidant and anaphylatoxin activity in vitro (Balint et al., Cancer Res. 44, 734, 1984).

20

The commercialisation of recombinant Protein A species as set forth in US6,399,750 which recombinant species is attached to the column matrix via a single thioester bond allowed for higher capacity protein A columns. As a concomittant disadvantage, the leakage rate of such recombinant Protein A matrices is often drastically increased in  
25 contrast to many traditional, multi-point attached natural Protein A matrices obtained by CNBr coupling. Protein A contaminant removal should therefore proceed without concomitant removal of complexed IgG.

Balint et al. (ibid.) demonstrated that such IgG-Protein complexes can be separated from  
30 uncomplexed IgG by gel filtration. Low through-put and loss in antibody yield are the disadvantages of this method.

US 4,983,722 teaches selective separation of contaminating protein A from a protein-A purified antibody preparation by absorbing the mixture to an anion exchanger material and to separate both components by sequentially eluting the antibodies and protein A under conditions of increasing ionic strength. This resolution method is highly dependent on the pI of the antibody which is specific and highly variable for a given antibody. Further, throughput is limited by the steepness of the salt gradient required for obtaining separation.

It is an object of the present invention to devise another method for separating protein A or protein A fragments from antibody, preferably an IgG, which method avoids the disadvantages of the prior art. According to the present invention, such object is solved according to the independent claims 1 and 9.

According to the present invention, a method of purifying an antibody is devised which method comprises the steps of: First, purifying an antibody by means of protein A affinity chromatography wherein the protein A is a native protein A or a functional derivative thereof.

Second, loading the purified antibody on an ion exchange material under conditions which allow for binding of the protein A or its functional derivative and third, collecting the antibody, preferably collecting at least 70%, more preferably collecting at least 80%, most preferably collecting at least 90% of the amount of antibody loaded onto the ion exchange material in the flow-through of the ion exchanger whilst any contaminant protein A or protein A derivative is bound to the ion exchange material.

Protein A is a cell surface protein found in *Staphylococcus aureus*. It has the property of binding the Fc region of a mammalian antibody, in particular of IgG class antibodies. Within a given class of antibodies, the affinity slightly varies with regard to species origin and antibody subclass or allotype (reviewed in Surolia, A. et al., 1982, Protein A: Nature's universal 'antibody', TIBS 7, 74-76; Langone et al., 1982, Protein A of *staphylococcus aureus* and related immunoglobulin receptors, *Advances in Immunology* 32:157-252). Protein A can be isolated directly from cultures of *S. aureus* that are secreting protein A or is more conveniently recombinantly expressed in *E.coli* (Lofdahl et al., 1983, *Proc. Natl. Acad. Sci. USA* 80:697-701). Its use for purification of antibodies, in particular monoclonal IgG, is amply described in the prior art (e.g. Langone et al., supra; Hjelm et

al, 1972; FEBS Lett. 28: 73-76). For use in protein A affinity chromatography, protein A is coupled to a solid matrix such as crosslinked, uncharged agarose (Sephacrose, freed from charged fraction of natural agarose), trisacryl, crosslinked dextrane or silica-based materials. Methods for such are commonly known in the art, e.g. coupling via primary amino functions of the protein to a CNBr-activated matrix. Protein A binds with high affinity and high specificity to the Fc portion of IgG, that is the C $\gamma$ 2-C $\gamma$ 3 interface region of IgG as described in Langone et al., 1982, supra. In particular, it binds strongly to the human allotypes or subclasses IgG1, IgG2, IgG3 and the mouse allotypes or subclasses IgG2a, IgG2b, IgG3. Protein A also exhibits an affinity for the Fab region of immunoglobulins that are encoded by the V<sub>H</sub> gene family, V<sub>H</sub> III (Sasso et al., 1991, J. Immunol, 61: 3026-3031; Hilson et al., J Exp. Med., 178: 331-336 (1993)). The sequence of the gene coding for protein A revealed two functionally distinct regions (Uhlen et al., J. Biol. Chem., 259: 1695-1702 (1984); Lofdahl et al., Proc. Natl. Acad. Sci.(USA), 80: 697-701 (1983)). The amino-terminal region contains five highly homologous IgG-binding domains (termed E, D, A, B and C), and the carboxy terminal region anchors the protein to the cell wall and membrane. All five IgG-binding domains of protein A bind to IgG via the Fc region, involving e.g. in human IgG-Fc residues 252-254, 433-435 and 311, as shown for the crystal structure in Deisenhofer et al. (1981, Biochemistry 20: 2361-2370) and in Sauer-Eriksson et al. (1995, Structure 3: 265-278) in case of the B-domain of protein A.

The finding of two essentially contiguous main binding sites in the Fc portion has been confirmed in the NMR-solution study of Gouda et al., 1998, Biochemistry 37: 129-136. In principle, each of the IgG-binding domains A to E of protein A is sufficient for binding to the Fc-portion of an IgG.

Further, certain alleles of the VH3 family in man have been found to mediate optionally binding of human Ig by protein A (Ibrahim et al., 1993, J. Immunol. 151:3597-3603; V-region mediated binding of human Ig by protein A). In the context of the present application, in another, separate object of the present invention, everything that has been said applying to Fc-region binding of antibody to protein A applies likewise to the binding of antibodies via such VH3 family protein A-binding allele in case that the Fc-region of such antibody did not allow on itself for high-affinity protein A binding. It may be considered an equivalent embodiment of the principal, Fc-based method of the present invention; the latter is further described in the subsequent sections.

An IgG antibody according to the present invention is to be understood as an antibody of such allotype that it can be bound to protein A in a high-affinity mode. Further, apart from the Fc portions of the antibody that are relevant for binding to protein A, such antibody must not correspond to a naturally occurring antibody. In particular in its variable chain regions portions, it can be an engineered chimeric or CDR-grafted antibody as are routinely devised in the art. An IgG-antibody according to the present invention is to be understood as an IgG-type antibody, in short.

A functional derivative of protein A according to the present invention is characterized by a binding constant of at least  $K=10^{-8}$  M, preferably  $K=10^{-9}$  M for the Fc portion of mouse IgG2a or human IgG1. An interaction compliant with such value for the binding constant is termed 'high affinity binding' in the present context. Preferably, such functional derivative of protein A comprises at least part of a functional IgG binding domain of wild-type protein A which domain is selected from the natural domains E,D,A,B, C or engineered mutants thereof which have retained IgG binding functionality. An example of such is the functional 59 aminoacid 'Z'-fragment of domain B of protein A which domain may be used for antibody purification as set forth in US 6013763. Preferably, however, an antibody binding fragment according to the present invention comprises at least two intact Fc binding domains as defined in this paragraph. An example of such are the recombinant protein A sequences disclosed e.g. in EP-282 308 and EP-284 368, both from Repligen Corporation.

Alone or in combination with a protein A or a functional protein A-fragment or derivative as defined in the preceding sections, further preferred are protein A fragments that are engineered to allow of single-point attachment. Single point attachment means that the protein moiety is attached via a single covalent bond to a chromatographic support material of the protein A affinity chromatography. Such single-point attachment by means of suitably reactive residues which further are ideally placed at an exposed amino acid position, namely in a loop, close to the N- or C-terminus or elsewhere on the outer circumference of the protein fold. Suitable reactive groups are e.g. sulfhydryl or amino functions. More preferably, such recombinant protein A or functional fragment thereof comprises a cysteine in its amino acid sequence. Most preferably, the cysteine is comprised in a segment that consists of the last 30 amino acids of the C-terminus of the amino acid

sequence of the recombinant protein A or functional fragment thereof. In a further preferred embodiment of such type, the recombinant protein A or functional fragment thereof is attached by at least 50% via a thioether sulphur bond to the chromatographic support or matrix material of the protein A-affinity chromatography medium. An example of such an embodiment is described e.g. in US 6399750 from Pharmacia and is commercially available under the brandnames of Streamline™ or MabSelect™ from Amersham-Biosciences, depending on the nature of the support matrix used. In the present context, thioether is to be understood narrowly as a –S– bonding scheme irrespective of chemical context, deviating in this regard from normal chemical language; it is possible, for instance, that said –S– ‘thioether’ bridge according to the present invention is part of a larger functional group such as e.g. a thioester or a mixed acetal, deviating in this regard in the context of the present application from the reactivity-based normal language of chemists. Preferably, the thioether bridge is a thioether bridge in its ordinary, narrow chemical meaning. Such bridging thioether group can be e.g. generated by reacting the 15 sulfhydryl-group of a cysteine residue of the protein A with an epoxide group harbored on the activated chromatographic support material. With a terminal cysteine residue, such reaction can be carried out under conditions suitable as to allow only for coupling of an exposed, unique sulfhydrylgroup of a protein as to result in single-point attachment of such protein only.

20

In a particularly preferred embodiment, the protein A or functional protein A derivative according to the present invention is the recombinant protein A disclosed in US 6399750 which comprises a juxta-terminal, engineered cysteine residue and is, preferably by at least 50%, more preferably by at least 70%, coupled to the chromatographic support material through the sulphur atom of said cysteine residue as the sole point of attachment. Further preferred, such coupling has been achieved by means of epoxide mediated activation, more preferably either by means of 1,4-bis-(2,3-epoxypropoxy) butane activation of e.g. an agarose matrix such as Sepharose Fast Flow (agarose beads crosslinked with epichlorohydrin, Amersham Biosciences, UK) or by means of epichlorohydrin activation of e.g. an agarose matrix such as Sepharose FF. Further preferred in combination with 30 afore said preferred embodiment according to this paragraph is that the first ion exchanger is an anion exchanger, in particular a quaternary amine-based anion exchanger such as Sepharose Q™ FF (Amersham-Biosciences/Pharmacia), most preferably it is an anion

exchanger having the functional exchanger group Q coupled to a matrix support which group Q is N,N,N-Trimethylamino-methyl, most preferably the anion exchanger is Sepharose Q<sup>TM</sup> FF from Pharmacia/Amersham Biosciences. The quarternary amino group is a strong exchanger which further is not susceptible to changes in pH of the loading/wash  
5 buffer. The fast flow exchanger matrix is based on 45-165 µm agarose beads having a high degree of crosslinking for higher physical stability; further sepharose is devoid of the charged, sulfated molecule fraction of natural agarose and does not allow for unspecific matrix adsorption of antibody, even under condition of high antibody loads. An example of such an embodiment can be found in the experimental section.

10

A contaminant protein A according to the present invention is any type of functional, IgG binding offspring of a protein A or a functional derivative thereof as defined above which is obtained upon eluting bound antibody from a protein A affinity chromatography column. Such contaminant protein A species may result e.g. from hydrolysis of peptide bonds  
15 which is very likely to occur by means of enzyme action in particular in industrial manufacturing. Protein A chromatography is applied as an early step in downstream processing when the crudely purified, fresh product solution still harbors considerable protease activity. Dying cells in the cell culture broth or cells disrupted in initial centrifugation or filtration steps are likely to have set free proteases; for regulatory  
20 purposes, supplementation of the cell culture broth with protease inhibitors prior or in the course of downstream processing is usually not accomplished, in contrast to biochemical research practice. Examples are Phenyl-methyl-sulfonyl-chloride (PMSF) or ε-caproic acid. Such chemical agents are undesirable as additives in the production of biopharmaceuticals. It is further possible that recombinant functional derivatives or  
25 fragments of protein A are less protease resistant than wild-type protein A, depending on the tertiary structure of the protein fold. Amino acid segments linking individual IgG binding domains might be exposed once the total number of binding domains is reduced. Interdomain contacts may possibly contribute to the stability of domain folding. It might also be that binding of antibody by protein A or said functional derivatives thereof  
30 influences or facilitates susceptibility to protease action, due to conformational changes induced upon binding of the antibody. Again, wild-type or full length protein A or functional, engineered fragments thereof might behave differently. Preferably, contaminant protein A according to the present invention still is functional, IgG binding protein and

thus is associated with the protein A-purified antibody when loaded onto the subsequent ion exchange separation medium according to the present invention. The high-affinity based association of contaminant protein A with the purified antibody is the reason why it is difficult to efficiently separate contaminant protein A from purified antibody.

5

Preferably, according to the present invention the antibody sought to be purified is harvested from a cell culture prior to purifying the antibody by means of protein A affinity chromatography. More preferably, said cell culture is a mammalian cell culture.

Mammalian cells have large compartments called lysosomes harboring degrading enzymes which are disrupted upon cell death or harvest. In particular, said cell culture may be a myeloma cell culture such as e.g. NS0 cells (Galfre, G. and Milstein, C. Methods Enzymology, 1981, 73,3). Myeloma cells are plasmacytoma cells, i.e. cells of lymphoid cell lineage. An exemplary NS0 cell line is e.g. cell line ECACC No. 85110503, freely available from the European Collection of Cell Cultures (ECACC), Centre for Applied Microbiology & Research, Salisbury, Wiltshire SP4 0JG, United Kingdom. NS0 have been found able to give rise to extremely high product yields, in particular if used for production of recombinant antibodies. In return, NS0 cells have been found to give reproducibly rise to much higher levels of contaminant protein A than other host cell types at least with certain protein A affinity chromatography systems employing recombinant, shortened fragments of wild-type protein A which recombinant protein A is possibly single-point attached protein A. An example of such is Streamline™ rProtein A affinity chromatography resin (Amersham Biosciences; essentially thioester single-point attached recombinant protein A as described in US 6,399,750). Levels of about or in excess of 1000 ng contaminant protein A/mg antibody could be obtained with Streamline™ rProtein A affinity columns. The method of the present invention distinguishes from the prior art in efficiently reducing contaminant protein A from such elevated levels to < 1 ng/mg antibody in a single, fast purification step, that is with a purification factor of about 1000x.

Further preferred is, alone or in combination with the preceding paragraph, that the antibody that is to be purified by means of protein A affinity chromatography is not treated as to inactivate proteases at or after harvest, more preferably is not in admixture with at least one exogenously supplemented protease inhibitor after harvest. Most preferably, said protease inhibitor is selected from the group consisting of PMSF, specific proteinase

30



inhibiting peptides as described in Laskowski et al., 1980, Protein inhibitors of proteinases, Ann. Rev. Biochem. 49, 593-626, and epsilon-caproic acid.

Operation of protein A affinity chromatography has been widely described in the technical literature and does not need to be further described. Another example apart from the above  
5 cited is e.g. Duhamel et al., J. Immunological Methods 31, (1979) 211-217, pH Gradient elution of human IgG1, IgG2 and IgG4 from protein A-Sepharose.

Preferably, the contaminant protein A is reduced to a concentration of < 10 ng/mg  
10 antibody, more preferably < 4 ng/mg antibody, most preferably < 1 ng/mg antibody in the flow-through of the first ion-exchanger, wherein antibody is preferably to be understood as to refer to IgG. The Elisa assay method for validation of these threshold values is described in detail in the experimental section; it should be noted that acidification of the sample to a  $\text{pH} \leq 4$ , preferably in the presence of a mild detergent, is crucial for accurate determination  
15 of the amount of leaked protein A. It goes without saying that this threshold is to be understood such as that the loading capacity of the first ion-exchanger for protein A binding is never exceeded, leading inevitably to break-through of contaminant protein A. A suitable Elisa-based method for assaying protein A or protein A fragments is described in US 4,983,722. Suitable anti-protein A antibodies are commercially available, e.g. from  
20 Sigma-Aldrich. In particular when using derivatives of protein A which derivatives have been engineered to harbor additional sulfhydryl groups, proper maintenance of the protein standard is important. It may be important to verify the monomeric character of such pure protein A derivative used as a standard for quantification of the test sample, since covalent di- or multimers formed via -S-S- bridges could lead to wrong results. Verification can be  
25 easily achieved by SDS-PAGE analysis under reducing and non-reducing conditions, as is customary in the art. Reduction of such protein A derivative- standard solution by means of DTT or beta-mercaptoethanol helps accordingly to circumvent errors of measurement in the ELISA-technique.

30 Further preferred, in the method according to the present invention at least 70%, more preferably at least 80%, most preferably at least 90% of the antibody loaded onto the first ion exchanger can be recovered in the flow-through of the ion-exchanger.

The first ion exchanger according to the present invention is an anion exchanger resin; protein A can be bound by both types of resin as described (EP-289 129 B1). The first ion exchanger or anion exchanger can be operated in the column mode at a certain flow rate or in batch operation mode, by submerging the ion exchange resin into the mildly agitated sample solution and further exchanging liquid media by filtration subsequently. According to the present invention and taking into account the pI of a given antibody, it is possible to define suitable conditions of pH and ionic strength for loading the first ion exchanger, which conditions result in retaining the antibody in the flow through whilst the protein A contaminant is bound and thus removed from the antibody. As has been said before, the method according to the present invention allows of faster separation of antibody from contaminant protein A. Examples of functional groups of such first, anion exchanger that are attached to a matrix support are e.g. primary, secondary, and particularly tertiary or quaternary amino groups such as aminoethyl, diethylaminoethyl, trimethylaminoethyl, trimethylaminomethyl and diethyl-(2-hydroxypropyl)-aminoethyl. Suitable chromatographic support matrixes for the anion exchanger are known in the art. Examples are agarose-based resins and beads, dextran beads, polystyrene beads and polystyrene/divinyl-benzene resins. Most preferably, the ion exchanger is a quaternary amine-based anion exchanger mounted on an agarose matrix such as e.g. Sepharose CL-6B or Sepharose Fast Flow (FF) from Amersham-Biosciences/Pharmacia. An example of such is Sepharose Q<sup>TM</sup> from Amersham-Biosciences/Pharmacia. Further preferred in conjunction with the use of a first anion exchanger is that the antibody according to the present invention is a monoclonal antibody that has an isoelectric point (pI) which is at least two pH units above, that is it is more basic than, the pI of the protein A used in the preceding protein A affinity chromatography step; e.g. whereas native protein A has a pI of about 5.0, Streamline recombinant protein A has a pI of about 4.5. Preferably, the antibody according to the present invention is a monoclonal antibody that has an isoelectric point (pI) which is at least 6.5 or above, more preferably is 7.0 or above, most preferably has an pI of at least 7.5 or above. It should be noted that this refers to the pI of the actually harvested and purified antibody, not the pI that can be simply predicted from the amino acid sequence alone. The actually purified antibody molecule may have undergone further modifications of the polypeptide backbone such as glycosylation, which modifications may add charged moieties and thus may have changed the pI of the molecule. Upon determination of pI for product antibody by means of isoelectric focusing (IEF), the

microheterogeneity of posttranslational processing of the antibody protein, e.g. a monoclonal antibody protein, leads to a wider pI-range for individual product antibody glycoprotein molecules, them resembling a smear in an IEF gel rather than a single band and thus a specific numeric value for at least the majority of product. According to the present invention, in such afore mentioned preferred embodiment, the 'pI of an antibody' refers to that share of antibody product molecules whose pI is within the preferred range of pI as specified above. All further definitions of this description, such as the %-proportion of antibody recovered after a given purification step, refer to said pI-compliant share of antibody only.

10

15

20

25

30

The mode of operation of a first anion exchanger according to the present invention requires buffer exchange of the acidic or neutralized eluate from the protein A affinity chromatography step with the equilibration buffer of the first anion exchanger. Equilibration buffer and loading buffer are identical in the method of the present invention. Commonly employed ultrafiltration devices such as sold by Amicon or Millipore can be expediently used for that purpose; those avoid the dilution effects whilst using e.g. low molecular weight porous filtration matrices such as Sephadex G-25. The equilibration buffer according to the present invention preferably has a salt concentration of a displacer salt such as e.g. sodium chloride in the range of 1 to 150 mM, more preferably of from 5 to 110 mM, most preferably of from 20 to 100 mM salt. The pH of the equilibration buffer is preferably in the range of pH 6.5 to pH 9.0, more preferably is in the range of pH 7.5 to pH 8.5, most preferably is in the range of pH 7.9 to pH 8.4. It should be kept in mind that N-terminal amino function of a protein has a pKs value of about 9.25, thus binding of contaminant protein A and any other already negatively charged protein to an anion exchanger will get stronger at more basic pH; for a given application, the pH of the loading buffer might need finetuning for optimal discrimination of binding and non-binding for a given pair of antibody and contaminant protein A having differing pI values and different content of cysteine and histidine side chains which may contribute to changes in charge within the selected ranges of pH. Further, a more basic pH interferes with proteinA-antibody interactions as will do any increase in ionic strength; likewise, ionic strength needs finetuning to balance prevention of binding of antibody with the need to abolish binding of contaminant

protein A. It goes without saying for the skilled artisan that the ionic strength of the buffer is usually inversely correlated with the pH value; the more strongly protein A gets bound to the anion exchanger depending on pH, the more salt is tolerated for preventing binding of antibody and for interfering with potential proteinA-antibody interactions. Thus, the above given ranges for pH and displacer salt thus are to be understood as to be correlated: The lower the pH, the less salt is found permissible within the above given preferred ranges for working the invention. Further salt freight is added by the pH buffering substance, further increasing the ionic strength of the solution. The ionic strength can be determined by measuring the conductivity of the equilibration buffer. The term 'conductivity' refers to the ability of an aqueous solution to conduct an electric current between two electrodes measures the total amount of ions further taking charge and ion motility into account. Therefore, with an increasing amount of ions present in the aqueous solution, the solution will have a higher conductivity. The unit of measurement for conductivity is mS/cm (milliSiemens/cm), and can be measured using a commercially available conductivity meter, e.g. from Topac Inc. (Hingham, MA/U.S.A.) or Honeywell. In the context of the present application, all numerical values pertain to the specific conductivity at 25°C. Preferably, the loading or equilibration buffer for the first anion exchange step has a conductivity 0.5-5 mS/cm, more preferably of from 1-3 mS/cm, most preferably of from 1.25-2.5 mS/cm. Ideally, it has a conductivity of about 2 mS/cm. Examples of suitable buffer salts can be found in Good, N.E. (1986, Biochemistry 5:467-476). E.g. Tris.HCl buffer or a sodium hydrogen phosphate buffer as customarily employed are suitable buffering substances. The concentration of the buffer substance is customarily in the range of e.g. 10-40 mM buffer salt. Amongst potential anion species useful for devising a buffer, those having lower specific strength of anion elution as compared to chloride, which property of low elution strength is approximately inversely correlated with the density of ionic charge and is approximately proportional to the ionic size, are preferred. Empirical comparisons of strength of anionic elution are tabulated in the standard textbooks of biochemistry. More preferably, the buffer substance according to the present invention is a phosphate buffer. Hydrogenphosphate has a low elution strength, in particular if employed at a pH at or below pH 8, and further excels by particularly low chaotropic properties.

Whereas batch mode operation is possible, column operation mode is preferred for the first anion exchanger step. In that case, a flow rate of about 10 to 60 ml/h can be advantageously employed. The loading concentration of antibody loaded can favorably be in the range of 10 to 30 mg antibody/ml exchange resin. It goes without saying that the use of extremely diluted samples would give rise to decreased yield of antibody, as is known to the skilled person. The antibody sought to be purified is collected in the low-through of the loading operation including about one column volume of wash with the same equilibration buffer. The pH of the flow-through may be adjusted to neutral pH for improving stability and preventing aggregation and/or precipitation of antibody protein.

After the first anion exchanger, the antibody is ready for use in applications or may be deemed to require further polishing by customary purification methods. In a further preferred embodiment, the first ion exchange step is followed by a second ion exchange step in which second step the antibody is loaded and bound by the second ion exchange medium and is eluted with a buffer other than the loading buffer, by means of increased salt and/or pH, as an essentially monomeric, non-aggregated antibody. 'Essentially' means less than 5% in this context. Preferably, alone or in combination with a preferred embodiment described in the preceding sections, the second ion exchanger is a cation exchanger. Such combination of a protein A chromatography step followed by a first anion exchanger and a second cation exchanger step is novel. It is well known that most trace contaminant proteins from cell culture broth have much lower pI values than antibodies, in particular IgG antibodies; cation exchange will therefore allow of efficient removal both of aggregated antibody and potential infectious agents such as virus capsids as well as of protein contaminants other than antibody. Due to speedy operation, highly efficient recovery of antibody after loading, binding to and elution from the co and high capacity of loading, it allows also of repeated, cyclic operation with a single batch of antibody with additive effect of the purification factor achieved in a single round of binding and elution. Preferably, the pH of the loading buffer is about pH 4 to 7, more preferably pH 4.01 to 6, most preferably pH 4.02 to 5.5. Further preferred, the the antibody is eluted from the cation exchanger with a salt gradient in the range of from 0.1 to 1.2 M salt, wherein the salt preferably is an alkaline metal salt, more preferably a lithium,

potassium or sodium salt. Preferably, elution takes place at a pH of from pH 7 to 8 in order to have maximum aggregate removal and minimal damage to antibody due to acidic conditions. Optionally preferred, the elution takes place at a an acid pH of from pH 4 to 7, more preferably 4.01 to 6 for maximizing removal of contaminant protein A; levels as low  
5 as < 0.4 ng/mg antibody can be realized in this way. This second cation exchanger step renders traditional gel filtration moot whilst allowing of high-capacity as well as fast operation as is typical for ion exchangers. Ion exchangers support loads of 10-30 mg antibody/ml resin. In a particularly preferred embodiment, the purification method of a first anion exchanger and a second cation exchanger step in the aftermath of protein A  
10 chromatography renders clinical grade antibody in the absence of a further, terminal size exclusion chromatography (SEC) step which SEC step would have a molecular weight cut off suitable for separating antibody aggregates and/or antibody-protein A complexes from monomeric antibody such as an normal IgG.

15 On a general note, the method of the present invention can not be exploited for antibodies that have been raised against protein A-borne epitopes. Such antibodies are disclaimed, though this is an obvious limitation to the skilled artisan.

The most appealing feature of the method of the present invention is that purifying  
20 antibody via an anion exchanger in a non-binding or flow-through mode, the capacity of the column is not all limiting the through-put of material; the capacity is only decisive with regard to minor amounts of contaminant protein A retain. This saves a lot of processing time and material resources whilst allowing for very efficient removal of protein A contaminant.

25

## Experiments

### 1. Protein A Elisa

Numerous Elisas for testing of protein A or recombinant protein A have been described  
5 (see US 4983722 and references described in there). For all work described below, a  
simple sandwich Elisa was used in which capture anti-protein A antibody coated on a flat-  
bottomed 96 well microtiter plate (Nunc™) retains the protein A. Bound protein A is then  
detected an a biotinylated anti-protein A detection antibody, which allows for further  
binding of streptavidin conjugated horseradish peroxidase (Amersham #RPN 1231).  
10 Commercially available anti-protein A rabbit antibody (raised against the natural *S. aureus*  
protein A) for capture is available from Sigma-Aldrich (#P-3775). It was this antibody  
which was used through-out this study. The detection rabbit antibody was equally  
purchased from Sigma-Aldrich (#3775). After coating the protein by unspecific adsorption  
process, the coated protein is used to retain protein A-specific protein A capture antibody  
15 which capture antibody is further detected with bioinylated rabbit anti-protein A and  
streptavidin-horseradish peroxidase. Tetramethyl benzidine is used as the chromogenic  
substrate. Samples of unknown concentration are read off against a standard curve using  
the very parent-protein A or -protein A derivative of the contaminant protein A sought to  
be detected. Coating at acidic pH as well as proper preparation of the standard has proven  
20 important. In particular for recombinant protein A's engineered to carry additional cysteine  
residue such as e.g. Streamline protein A™ (Amersham Biosciences, formerly Pharmacia),  
the standard solution was found to require pretreatment with a reducing sulfhydryl agent to  
ensure monomeric state of the protein standard solution.

Wild-type protein A standard, in contrast, is commercially available from a number of  
25 companies, e.g. Sigma-Aldrich/Switzerland (#P6031) or Pharmacia (#17-0770-01) and  
does not require such pretreatment. For the below described experiments observing leakage  
of contaminant protein A from Streamline™ matrix, samples of unconjugated recombinant  
protein A obtained from the manufacturer were used as a standard.

### 1.1 Pretreatment of Cys-enriched protein A-standard

Pure recombinant protein A-Cys as commercially in the Streamline™ protein A affinity chromatography (Amersham Biosciences) column material was obtained freeze-dried from Pharmacia/Amersham Biosciences. Up to 20 mg/ml protein were dissolved in 0.1M Tris pH 8 containing 0.5 M NaCl, 1mM EDTA and 20 mM dithioerythritol, incubated for 15-30 min. at room temperature and desalted with a disposable PD-10 gel filtration column (Amersham Biosciences). All buffers used for handling the standard solution before coating should be N<sub>2</sub>-treated to prevent oxidation of the thiol groups. Preparation of the protein standard was carried out at best immediately prior to use of the standard for coating the microtiter plates. Optionally, a 1mg/ml stock solution was prepared and kept at -65°C in a freezer; after thawing, monomeric character of protein A was assayed from an aliquot loaded on non-reducing SDS-PAGE. The concentration of protein standard was determined by Bradford assay (Bradford et al., 1976, Anal. Biochem. 72:248-254; Splittgerber et al., 1989, Anal. Biochem. 179:198-201) as well as by automated amino acid analysis. The result of such pretreatment is shown in Fig. 1 by means of non-reducing 10% SDS-PAGE for a staphylococcal protein A standard (lane 1: native protein A; lane 2: after pretreatment) and pure, uncoupled Streamline™ recombinant protein A (provided by courtesy of Pharmacia, now Amersham-Biosciences; lane 4: native recombinant protein A; lane 5: after pretreatment). Lane 1 is a molecular weight marker with the corresponding molecular masses being denoted on the vertical axis. The recombinant protein A from Pharmacia harboring an additional Cys residue shifts after reduction to lower molecular weight; a monomeric band at about 34 kD is preserved and much more intense, stemming obviously from dissociation of disulfide bridged dimers.

### 1.2 Elisa

#### 1.2.1 Preparation of sample

By two dilution steps, 1: 200.000 dilution of the 1 mg/ml protein A standard stock solution was prepared to provide the top standard at 50 ng/ml. Thereof, dilutions down to 0.2 ng/ml were prepared for assaying the standard curve. Further, the dilutions of the standard ('spinking solutions') were used for spiking of duplicate



product samples to be tested in order to exclude presence of interfering substances in the sample.

For final product sample testing, every sample is divided into 2 equal volumes of 500  $\mu$ l. One is spiked with the 1000 ng/ml spiking solution, or the 10  $\mu$ g/ml solution if appropriate, to give a final protein A content of 10 ng protein A per mg of antibody. The other half is spiked with the same volume of sample buffer; thus the dilution factor of the product sample due to spiking is accounted for. Both types of preparation will be referred to as 'spiked sample' in the following. The sample buffer was made up from 7.51 g Glycine (base), 5.84 g NaCl, 0.5 ml Triton X-100 to a volume of 1 L with deionized or bidestillated water.

For optimal accuracy measurements, the antibody concentrations in the samples were pre-determined by customary Elisa's well known in the art. A further standard solution was spiked with an equal amount of a known standard antibody of comparable constant region affinity for protein A, to determine efficiency of the acidification step and to unravel any potential systematic error introduced by antibody binding to and thus scavenging protein A from capture in the assay.

**Acidification:** To 450  $\mu$ l of spiked sample or standard is added 200  $\mu$ l of 0.2 M citrate/0.05% Triton X-100 buffer at pH 3.0. All samples were done in triplicate. Further, dilutions of sample were prepared and tested in triplicate since the assay works optimal for antibody concentrations being in the range of 1 mg/ml and 0.2 mg/ml. The acidification step is crucial in the present assay to liberate contaminant protein A or A fragments which were otherwise bound to the excess of antibody present in the sample solution.

### 1.2.2 Coating of microtiter plates with antibody

Coating buffer was made up from 1.59 g/L Na<sub>2</sub>CO<sub>3</sub>, 2.93 g/L NaHCO<sub>3</sub> and 0.20 g/L sodium azide. The pH of the buffer was adjusted to pH 9.6. Add 100  $\mu$ l antibody solution per well comprising antibody in an amount sufficient as not to show saturation for the protein A standard. Cover plate with plastic film and place in humidity chamber. Incubate at 37°C overnight for approximately 18 hours. Rinse all wells 3 times with at least 300  $\mu$ l washing buffer [NaCl 5.8 g/L, Na<sub>2</sub>HPO<sub>4</sub> 1.15

g/L,  $\text{NaH}_2\text{PO}_4 \cdot \text{H}_2\text{O}$  0.26 g/L, EDTA 3.7 g/L, Tween-20 0.2 g/L, butanol 10 ml/L, pH 7.2], and tap dry. Add 250  $\mu\text{l}$  blocking buffer [coating buffer with 0.5 % casein hammarsten] to each well and incubate for 2 hours at ambient temperature on a benchtop orbital shaker (speed 120 rpm). Rinse all wells three times with at least 300  $\mu\text{l}$  washing buffer, and tap dry.

### 1.2.3 Incubation of sample and detection

Plate out standards and samples including any spiked samples with 100  $\mu\text{l}$  /well.

Cover plate with plastic film and incubate for 90 minutes at ambient temperature on an orbital benchtop shaker. Rinse all wells three times with at least 300  $\mu\text{l}$  washing buffer, and tap dry. Dilute biotinylated rabbit anti-protein A at the previously determined optimal dilution. Add 100  $\mu\text{l}$ /well, cover plate with plastic film and incubate for 90 minutes at ambient temperature on an orbital shaker. Repeat rinsing.

Dilute streptavidin-horseradish peroxidase at the previously-determined optimal dilution using conjugate buffer [ $\text{Na}_2\text{HPO}_4$  1.15 g/L,  $\text{NaCl}$  5.84 g/L,  $\text{NaH}_2\text{PO}_4 \cdot \text{H}_2\text{O}$  0.26 g/L, EDTA 3.73 g/L, Triton X-100 0.05% (v/v), pH 7.2]. Add 100  $\mu\text{l}$ /well, cover plate in plastic film and incubate for 45 minutes at ambient temperatur on an orbital shaker. Repeat rinsing. Add 100  $\mu\text{l}$  freshly-prepared tetramethyl-benzidine (TMB, ICN product number #980502) substrate solution. The substrate solution is prepared like this: A stock solution is prepared by dissolving 10 mg TMB in 1 ml DMSO. 10  $\mu\text{l}$  of that stock, further 10  $\mu\text{l}$  of  $\text{H}_2\text{O}_2$  are added to a 2.05 % (w/w) sodium acetate aqueous solution that was adjusted to pH 6.0 with 0.5 M citric acid. It goes without saying that all water used for preparing any reagent of the assay is of highest quality, that is deionized ultrapure or at least bidestillated water.

The substrate solution is incubated at ambient temperature for 8-11 minutes on a shaker. The reaction is then stopped by adding 50  $\mu\text{l}$  per well of stopping solution [13%  $\text{H}_2\text{SO}_4$ ]. Within 10 min. after addition of the stopping solution, the absorbance of the wells at wavelength 450 nm is determined on a plate-reading spectrophotometer.

The detection limit for such Elisa is 0.2 ng/ml Protein A, with a working range of

from 0.2 to 50 ng/ml. The interassay variability is less than 10%.

Fig. 2 shows the levels of leaked recombinant protein A in antibody eluates from Streamline™ recombinant protein A chromatography with single-point attached protein A in thioether linkage. The cycle number refers to repeated use after elution with 1 M sodium chloride and re-equilibration. Whereas leakage from cell culture broth from hybridoma cell culture was typically in the order of 500 ppm, other cell types gave levels as high as 1000 ppm. An overview on the rate of leakage from differently sourced matrices is given in Table 1; chromatography was performed according to manufacturer' instruction.

**Table 1**

Matrix	Supplier	Coupling chemistry	Typical leakage p.p.m	Working Capacity (mgml <sup>-1</sup> )	Flow Rate (cmh <sup>-1</sup> )
Native Protein A Sepharose 4FF	Amersham-Biosciences	Multi-point attached CNBr	10 -20	5-20	30 - 300
rmp Protein A Sepharose	Amersham-Biosciences	Multi point attached	10 -20	5-20	30 - 300
Poros A High Capacity	Applied Biosystems	Multi point attached	10 - 50	10	500 - 1000
Protein A Ceramic HyperD	Biosepra	Multi point attached	Up to 300	10-20	200 - 500
rProtein A Sepharose	Amersham-Biosciences	Single point attached Thioether linkage	50 -1000	20-40	30 - 300
MabSelect	Amersham-Biosciences	Single point attached Thioether linkage	50 -1000	20-40	500
STREAMLINE rProtein A	Amersham-Biosciences	Single point attached Thioether linkage	50 -1000	20-40	200 - 400

Fig. 3 further provides data on insubstantially reduced leakage of contaminant protein A during repeated runs of the protein A affinity chromatography with the same affinity matrix material ; wild-type protein A multipoint-attached Sepharose 4 FF (Amersham-Biosciences) was repeatedly used as described in section 2.1 below and the level of contaminant protein A in the eluate, before any further processing

of eluate, was determined by Elisa as described above.

## 2. Protein A and Sepharose Q chromatography

### 5 2.1 Protein A affinity chromatography with Streamline™

Cell culture supernatant from a NS0 myeloma cell culture was crudely purified by centrifugation and in depth filtration and concentrated by ultrafiltration; ultrafiltration was also used to exchange the culture fluid to PBS pH 7.5. The titer of the antibody #5 produced by the cells was 0.2 mg/ml, a total of 1 L buffer-exchanged supernatant was loaded. The pI of the monoclonal antibody #5 was 8.5. 10 The protein A Streamline™ column (5.0 ml volume) was previously equilibrated with 10 column volumes of 50 mM glycine/glycinate pH 8.8, 4.3 M NaCl; flow rate was at 200 cm/h. For loading, the column was operated at a flow rate of 50 cmh<sup>-1</sup>; loading capacity was about 20 mg/ml matrix material). Before elution, the 15 column was washed with at least 10 column volumes of glycine equilibration buffer supplemented with additional 200 mM NaCl and 0.1% Tween-20. Elution was achieved with elution buffer made up of 0.1 M glycine/HCl pH 4.0 buffer. Immediately after elution, fractions of eluate comprising the antibody peak were neutralized with an adequate aliquot of 0.5 M TrisHCl pH 7.5 and buffer 20 exchanged with an Amicon diafiltration device with loading/equilibration buffer (10mM Tris/HCl pH 8.0, 50 mM NaCl) of the present invention for the subsequent anion exchanger step for preventing longer exposure to acidic pH. The antibody concentration and the contaminant protein A concentration were determined as described above. The level of contaminant protein A in the eluate 25 amounted to 1434 ng/mg antibody before and amounted to 1650 ng/mg antibody after diafiltration. The recovery of antibody based on the titer of the buffer exchanged supernatant solution prior to loading was 81%; the concentration of antibody in the diafiltrated solution was 3.6 mg/ml.

## 2.2 Q-Sepharose FF anion exchange step in non-binding mode

The purified antibody from section 2.1 was further processed as described: A 5.0 ml Q-Sepharose FF column (Amersham-Biosciences) was packed 10 ml of 0.1 M NaOH, followed by 2 column volumes of 0.1 M Tris pH 8, and equilibrated in 10 column volumes of 10 mM Tris pH 8/50 mM NaCl, at a flow rate of 75 cm/h. After equilibration, the flow rate was reduced to 50 cm/h. 6 ml of the diafiltrated antibody solution was loaded onto the column and the flow-through was collected for further processing; the flow-through was continued to be collected until, after having loaded the column with the initial 6 ml and having continued thereafter with pure loading or equilibration buffer 10 mM Tris pH 8, 50 mM NaCl, the absorption of the flow-through monitored at 280 nm was back to baseline. The total recovery of antibody in the flow-through was 23 mg antibody (87%). Determination of the level of contaminant protein A resulted in < 3 ng/mg antibody.

Further processing of this Q-Sepharose purified antibody batch by gel filtration (size exclusion chromatography, SEC) over Sephacryl S-300 in 10 mM Phosphate pH 7.0, 140 mM NaCl buffer at a flow rate of 10 cm/h with a loading ratio of 15 mg antibody per ml gel was found not change this trace contaminant protein A level substantially any more. By experience, SEC may be used to further reduce levels of about 30-100 ng/mg contaminant protein A to about 1-5 ng/mg. Thus SEC has a very low purification factor with regard to trace amounts of protein A, possibly accounting for affinity interactions in between antibody and contaminant A. However, due to the unavoidable dilution of sample and slow processing with allows for some decay of the antibody protein, SEC will allow for 70% recovery only of the amount of antibody loaded. This means SEC will unavoidably result in loss of material whilst requiring much time.

The Q-Sepharose column was recycled for further use by separate elution in 2M NaCl and further equilibration as described above.

### 3. Protein A and Sepharose Q purification with subsequent cation exchange step

In a further experiment, the antibody from exp. 2.2. purified in a non-binding mode by Q-Sepharose anion exchange was used. Instead of testing a further, final SEC purification step, the antibody harvested in the flow-through of the Sepharose Q column was subjected to a second cation exchanger step with a SP-Sepharose FF (SP=Sulphopropyl-) matrix from Amersham-Biosciences. The SP-Sepharose FF allowed of a flow rate of 100 cm/h with a reproducible yield of 93% antibody after loading, washing and elution of the antibody from the cation exchanger.

For loading, the pH of the antibody solution obtained after Sepharose Q purification step was adjusted to pH 4.5-5.0 with 50 mM acetate buffer pH 4.5. The loading capacity was set with 10 mg/ml matrix material at a conductivity of load of 17 mS/cm. The 50 mM acetate buffer was further used for washing to baseline. A 50 mM Na acetate pH 4.5, 1 M NaCl high salt buffer was used for elution of antibody; monomeric antibody eluted first, whereas aggregates used to elute in the tail fractions at high ionic strength. Use of a less steep salt gradient by implementation of a salt gradient in the elution buffer before pumping on the column is equally feasible; direct application of a high-salt buffer results in less diluted antibody and consequently more precisely sampling and shorter times of residence in the acidic solution. After elution, the acidic buffer was quickly exchanged for PBS pH 7.5. The level of contaminant protein A in the pooled eluate was determined with < 0.4 ng/mg antibody, the antibody was shown to be >99% monomeric by means of size exclusion HPLC.

#### **4. Streamline™ Protein A affinity chromatography with custom-made, multipoint -attached protein A**

5 This multipoint-attached Streamline™ protein A-affinity matrix was custom made and supplied by Pharmacia Biotech (now Amersham-Pharmacia). It was made up by the manufacturer by coupling the same 34 kD Streamline™ -type recombinant protein A having a terminal Cys residue to the same Sepharose matrix material, but used traditional CNBr chemistry for activation and coupling instead of epoxide-  
10 mediated activation and selective reaction conditions for coupling of -SH groups only (see product information from manufacturer). The method of exp. 2.1 was repeated and the level of contaminant protein A was determined with 353 ng/mg antibody. This means that the mode of coupling partly accounts for increased protein leakage from high-capacity, single-point attached recombinant protein A  
15 affinity matrices; the modifications in amino acid sequence introduced into such recombinant protein A as compared to full-length wild-type protein A contribute considerably to increased protein leakage, too.

#### **5. Parallel comparison of methods: Comparison with Miles Method (US4,983,722)**

20 The Miles Patent (No: 4,983,722) claims that binding DEAE Sepharose used as a second chromatography step with a salt gradient (0.025M to 0.25M NaCl) for elution can reduce the leached Protein A content in the eluate to less than 15ng/mg  
25 antibody (range of protein A was 0.9 to 14 ng/mg of antibody).

30

35

**Table 2:**

**Comparison of Protein A residues in eluate samples of 6A1 Antibody purified on single and multipoint attached Protein A affinity matrices**

Matrix	Sample	Protein A levels (ng/mg)
rProtein A Sepharose (single point attached)	Protein A eluate	20.2
rpm Protein A Sepharose (multi-point attached)	Protein A eluate	2.16
Native Protein A Sepharose (multipoint attached)	Protein A eluate	<2.0

The aim of these experiments was to confirm these results using MabSelect ( new single point attached rProtein A matrix) with a lower pI antibody (pI 6.5-7.5), and to directly compare the non-binding Q-Sepharose method (using different equilibration/loading buffers) with the Miles Patent method.

**Method applied:**

The purification of 6A1 antibody (pI 6.5 – 7.5) included two chromatography steps consisting of MabSelect Protein A step followed by Q-Sepharose anion exchange chromatography (non-binding), or DEAE Sepharose chromatography (binding) step. See L0 9007 and L0 9375

**MabSelect Protein A Chromatography:**

Column matrix Mab Select recombinant Protein A (single point attached rPA)  
Column dimensions 1.6 cm internal diameter x 15 cm bed height



	Column volume	30 mL
	Operational flow rate	500 cm/hr (16.80mL/min)
	Clean	6M guanidine HCL (2 column volumes)
	Loading capacity	35 mg /ml matrix
5	Equilibration	50mM glycine/glycinate pH 8.0/250mM NaCL (8 column volumes)
	Post load wash	50mM glycine/glycinate pH 8.0/250mM NaCL (8 column volumes)
	Elution buffer	100mM glycine pH 3.50 (6 column volumes)
10	Wash	100mM Citric acid pH 2.1 (2 column volumes)

The culture supernatant containing 6A1 antibody was purified on a MabSelect column (30ml), connected to an AKTA FPLC system. The conditions used were as described in the table above. The antibody was eluted using 0.1M glycine pH 3.5. Following elution  
 15 the eluate pH was adjusted to pH 7.0, and then the eluate sample was divided into 5 aliquots; each aliquot was then diafiltered into a different buffer for anion exchange chromatography.

The first aliquot was diafiltered into 50mMTrisHCl pH8 /75mMNaCl for Q-Sepharose chromatography run 1. The second aliquot was diafiltered into 50mMTrisHCl pH8 /100  
 20 mMNaCl for Q-Sepharose chromatography Run 2. The third aliquot was diafiltered into 20mM sodium phosphate pH6.5 /80 mM NaCl for Q-Sepharose chromatography Run 3. Aliquots four and five were buffer exchanged into 25mMTris HCl pH 8.0/25 mMNaCl for evaluation of binding DEAE Sepharose method described in Miles patent. The difference  
 25 between Runs 4 & 5 is that in Run 4 the main peak was collected as one fraction and diafiltered into standard phosphate buffered saline prior to analysis whereas in Run 5, the elution peak was fractionated and dialysed into a phosphate buffer prepared as described in the Miles Patent.

The conditions for each of the five column runs are described below:

### 30 **Q- Sepharose Chromatography: Run 1**

	Column matrix	Q-Sepharose Fast Flow
	Column dimensions	1.6 cm internal diameter x 8 cm bed height
	Column volume	16 mL
35	Column preparation	Packed in 0.1 M Sodium Hydroxide at 150 cm/hr
	Operational flow rate	100 cm/hr (3.35mL/min)
	Clean	0.1M Sodium Hydroxide (2 column volumes)
	Loading capacity	15 mg /ml matrix
	Equilibration	50mM TrisHCl pH 8.0/75mM NaCl (8 column volumes)
40	Post load wash	50mM TrisHCl pH 8.0/75mM NaCl (5 column volumes)
	Strip buffer	2 M Sodium Chloride (2 column volumes)

Wash 0.1M Sodium Hydroxide (2 column volumes)

### **Q- Sepharose Chromatography: Run 2**

- |    |                       |   |
|----|-----------------------|---|
| 5  | Column matrix         | Q-Sepharose Fast Flow                             |
|    | Column dimensions     | 1.6 cm internal diameter x 8 cm bed height        |
|    | Column volume         | 16 mL   |
|    | Column preparation    | Packed in 0.1 M Sodium Hydroxide at 150 cm/hr     |
|    | Operational flow rate | 100 cm/hr (3.35mL/min)                            |
| 10 | Clean                 | 0.1M Sodium Hydroxide (2 column volumes)          |
|    | Loading capacity      | 7.5 mg /ml matrix                                 |
|    | Equilibration         | 50mM TrisHCl pH 8.0/100mM NaCl (8 column volumes) |
|    | Post load wash        | 50mM TrisHCl pH 8.0/100mM NaCl (5 column volumes) |
|    | Strip buffer          | 2 M Sodium Chloride (2 column volumes)            |
| 15 | Wash                  | 0.1M Sodium Hydroxide (2 column volumes)          |

### **Q- Sepharose Chromatography: Run 3**

- |    |                       |  |
|----|-----------------------|--|
|    | Column matrix         | Q-Sepharose Fast Flow  |
| 20 | Column dimensions     | 1.6 cm internal diameter x 8 cm bed height                   |
|    | Column volume         | 16 mL  |
|    | Column preparation    | Packed in 0.1 M Sodium Hydroxide at 150 cm/hr                |
|    | Operational flow rate | 100 cm/hr (3.35mL/min)                                       |
|    | Clean                 | 0.1M Sodium Hydroxide (2 column volumes)                     |
| 25 | Loading capacity      | 7.5 mg /ml matrix  |
|    | Equilibration         | 20mM Sodium phosphate pH 6.5/80mM NaCl                       |
|    | Post load wash        | 20mM Sodium phosphate pH 6.5/80mM NaCl<br>(5 column volumes) |
|    | Strip buffer          | 2M Sodium Chloride (2 column volumes)                        |
| 30 | Wash                  | 0.1M Sodium Hydroxide (2 column volumes)                     |

### **DEAE Sepharose: Run 4**

- |    |                       |   |
|----|-----------------------|---|
|    | Column matrix         | DEAE Sepharose  |
| 35 | Column dimensions     | 1.6 cm internal diameter x 8 cm bed height  |
|    | Column volume         | 16 mL   |
|    | Column preparation    | Packed in equilibration buffer at 150 cm/hr   |
|    | Operational flow rate | 100 cm/hr (3.35mL/min)  |
|    | Clean                 | 0.1M Sodium Hydroxide (2 column volumes)  |
| 40 | Loading capacity      | 7.5 mg /ml matrix   |
|    | Equilibration         | 25mM TrisHCl pH 8.6/25mM NaCl<br>(8 column volumes)                                       |
|    | Post load wash        | 25mM TrisHCl pH 8.6/25mM NaCl<br>(5 column volumes)                                       |
| 45 | Elution buffer        | 25mM TrisHCl pH 8.6/25mM NaCl To 25mM TrisHCl pH<br>8.6/250mM NaCl<br>(10 column volumes) |
|    | Wash                  | 2M Sodium Chloride<br>(2 column volumes)  |

**DEAE Sepharose binding method: Run 5 (Miles method)**

5	Column matrix	DEAE Sepharose
	Column dimensions	1.6 cm internal diameter x 8 cm bed height
	Column volume	16 mL
	Column preparation	Packed in equilibration buffer at 150 cm/hr
10	Operational flow rate	100 cm/hr (3.35mL/min)
	Clean	0.1M Sodium Hydroxide (2 column volumes)
	Loading capacity	7.5 mg /ml matrix
	Equilibration	25mM TrisHCl pH 8.6/25mM NaCl (8 column volumes)
15	Post load wash	25mM TrisHCl pH 8.6/25mM NaCl (5 column volumes)
	Elution buffer	25mM TrisHCl pH 8.6/25mM NaCl To 25mM TrisHCl pH 8.6/250mM NaCl (10 column volumes)
20	Wash	2M Sodium Chloride (2 column volumes)

The properties of the different buffers used in this study are shown in Table 3.

The elution profiles for each of the anion exchange chromatography runs are shown in

25 Figures 2 – 5.

Eluate samples generated from the 5 ion exchange runs were assayed for Protein A levels in the rPA ELISA. The results are shown in Table 4.

**Table 3:**

30 **Buffers used in this study**

Equilibration Buffer	Run number	Conductivity (ms/cm)	Resin	pH
50mM TrisHCl pH 8.0 / 75mM NaCl	1	10.74	Q-Sepharose (non-binding)	8.00
50mM TrisHCl pH 8.0 / 100mM NaCl	2	12.85	Q-Sepharose (non-binding)	8.01
20mM Sodium phosphate pH 6.5 / 80mM NaCl	3	10.20	Q-Sepharose (non-binding)	6.50
25mM TrisHCl pH 8.6/25mM NaCl	4/5	3.35	DEAE-Sepharose (binding)	8.60
25mM TrisHCl pH 8.6/250 mM NaCl*	4/5	24.54	DEAE-Sepharose (binding)	8.61

\*Gradient elution buffer

5 Fractions across the elution profile of DEAE-Sepharose Run 5 (Miles method) were collected and analysed in the rProtein A ELISA; the results are shown in Table 5.

**Table 4: rProtein A ELISA Results:**

\*Where CV's denotes column volumes

Sample ID	rProtein A levels (ng/mg)	Antibody concentration (mg /ml)	% Recovery	Elution Volumes (CV's)*
Q-Sepharose eluate <b>Run 1</b>	< 0.4	1.42	82	4.5
Q-Sepharose eluate <b>Run 2</b>	2.94	1.49	70	3.5
Q-Sepharose eluate <b>Run 3</b>	0.73	1.86	85	3.4
DEAE Sepharose eluate pool <b>Run 4</b> (pool of all fractions)	1.72	2.16	75	2.5
DEAE Sepharose eluate pool (Miles Method) <b>Run 5</b> (pool of fractions 2 to 6)	1.55	1.83	73	3

10

**Table 5: Levels of rProtein A in Eluate fractions across the elution peak obtained during binding DEAE-Sepharose separation (Miles Method); Run 5.**

15


Fraction Number	rProtein A levels (ng/mg)	Absorbance (A <sub>280</sub> )
1	3.33	0.018
2	0.4	0.108
3	0.4	0.22
4	0.4	0.169
5	2.01	0.092
6	16.7	0.042
7	6.38	0.016

The highest recovery (85%) and best clearance of rProtein A for this antibody (6A1; pI 6.5 – 7.5) was obtained under non-binding conditions on Q-Sepharose using 20 mM sodium phosphate pH 6.5 / 80 mM NaCl buffer (corresponding to Run 3). Run 1 also showed good recovery (82%) and rProtein A clearance however, the elution volume for this run was significantly higher than expected for a non-binding method; suggesting partial retardation of the antibody on the column in this buffer system. Increasing the NaCl concentration (Run 2) resulted in lower rProtein A clearance, hence the buffer system used in Run 3 was more appropriate for this antibody. It has been our previous observation that the buffer system used in Run 1 is more appropriate for high pI antibodies and that used in Run 3 is particularly useful for neutral or slightly acidic antibodies. These experiments were done at similar capacities (7.5 mg/ml resin) we would expect to be able to use this non-binding method at much higher capacities (>30 mg/ml). We would expect this non-binding method to be applicable to many anion exchangers for example Q-Hyper D in addition to anion exchange membrane adsorbers (such as Mustang Q, Intercept Q and Sartobind Q). We would also expect this process to be more applicable to large scale production compared to the Miles method as higher capacities etc can be applied.

In the case of Run 5 (Miles method), fractionation of rProtein A was observed across the main elution peak as shown in table 5 . Careful pooling of fractions is therefore required to ensure good clearance of rProtein A. This had an impact on recovery (73%) and even in this case did not give as good clearance as obtained with the non-binding method. For the Miles method therefore it is more difficult to achieve good clearance and high recovery for cell lines/antibodies in cases where very high leakage is observed (such as that commonly obtained with single point attached matrices)

Data from **Run 5** is representative of the conditions described in the Miles patent.

An overview of method comparisons and the data obtained is shown in Table 6, below.

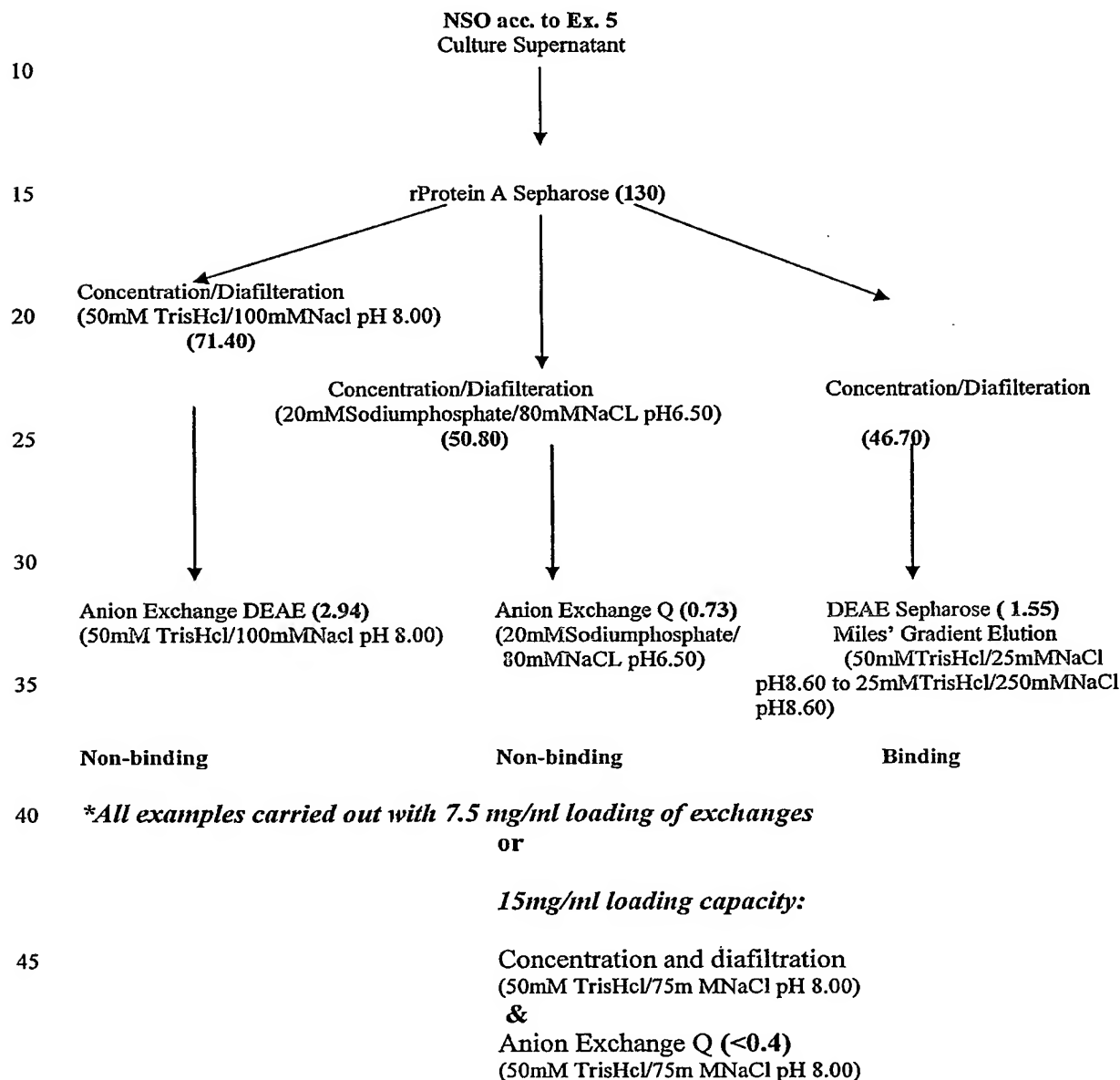**Table 6****Summary of rProtein A Levels at Different Stages of Antibody Purification.**

**Note:** Levels of rProtein A are shown in brackets ng/mg; note that not all NSO clonal cell lines supernatants give similar contamination levels of protein A.

## 6. Purification of a high pI Antibody

A high pI antibody (pI 9.0-9.3) was purified using Protein A Affinity Chromatography (MabSelect – single point attached recombinant Protein A matrix), followed by Q-

- 5 Sepharose anion exchange chromatography (under non binding conditions; for removal of trace contaminants) followed by SP-Sepharose cation exchange chromatography (under binding conditions for removal of aggregates).


10

15

20

25

30

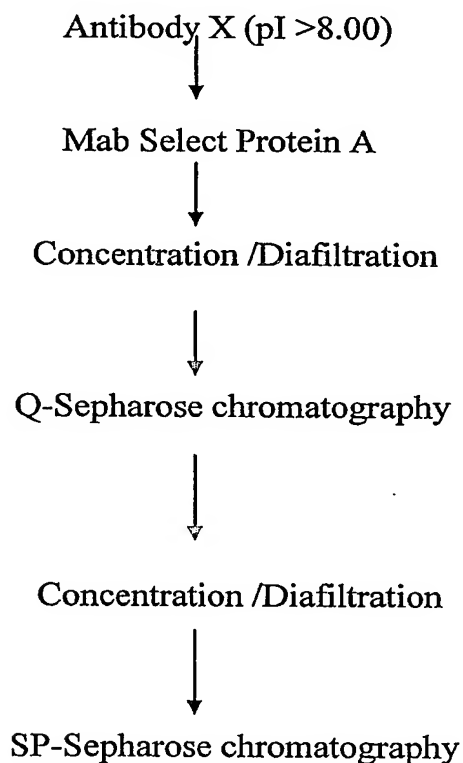

## Experimental Materials and Methods

35

### MabSelect Protein A Chromatography:

Column matrix	Mab Select recombinant Protein A (single point attached rPA)
Column dimensions	1.6 cm internal diameter x 15 cm bed height
40 Column volume	30 mL
Operational flow rate	500 cm/hr (16.80mL/min)
Clean	6M guanidine HCL (2 column volumes)

- 5      Loading capacity      35 mg /ml matrix  
      Equilibration      50mM glycine/glycinate pH 8.0/250mM NaCL  
                                   (8 column volumes)  
      Post load wash      50mM glycine/glycinate pH 8.0/250mM NaCL  
                                   (8 column volumes)  
      Elution buffer      100mM glycine pH 3.50 (6 column volumes)  
      Wash                  100mM Citric acid pH 2.1 (2 column volumes)
- 10      Culture supernatant containing high pI antibody was purified on a MabSelect Protein A Affinity column (30ml), connected to an AKTA FPLC system. The conditions used were as described in the table above. The antibody was eluted using 0.1M glycine pH 3.5. Following elution, the eluate was held at pH 3.69 (no adjustment required) for 60 min (low pH virus inactivation step), and then neutralised to pH 8 using 2 M Tris Base. Three cycles
- 15      on Protein A were performed; product recovery was determined by A280nm and is shown in Table 7 for each cycle.

**Table 7: % Recovery on Mab Select Protein A Column**

Cycle Number	% Recovery
1	81
2	81
3	80

- 20      After MabSelect Protein A chromatography the eluates from each of the three cycles were pooled together and buffer exchanged into 25mM Tris pH 8.0 (Q-Sepharose equilibration buffer) using an Amicon stirred cell concentrator fitted with 10kDa Millipore membrane.

25      **Q- Sepharose Chromatography:**

- Column matrix      Q-Sepharose Fast Flow  
      Column dimensions      1.6 cm internal diameter x 15 cm bed height  
      Column volume      30 mL  
      30      Column preparation      Packed in 0.1 M Sodium Hydroxide at 225 cm/hr  
      Operational flow rate      150 cm/hr (5.0 mL/min)  
      Clean                  0.1M Sodium Hydroxide (2 column volumes)  
      Loading capacity      40 mg /ml matrix  
      Equilibration      20mM Tris pH 8.0 (8 column volumes)  
      35      Post load wash      20mM Tris pH 8.0 (5 column volumes)  
      Strip buffer      20mM Tris pH 8.0 /2M NaCl (2 column volumes)  
      Wash                  0.1M Sodium Hydroxide (2 column volumes)



40ml of concentrated/diafiltered MabSelect Protein A eluate was loaded on Q-Sepharose column at a loading capacity of 40mg/ml matrix. The column was operated in a non-binding mode and the unbound fraction containing the antibody was collected. The recovery on this step was 69% by A280. This is slightly lower than obtained under these conditions for this antibody and may be due to inaccurate estimation of the load volume due to hold up volume in the FPLC sample pump.

Following Q-Sepharose chromatography, the unbound fraction was concentrated to 13.98 mg/ml and diafiltered into SP-Sepharose equilibration buffer (25mM Sodium acetate pH5.0/25mM NaCl) using an Amicon stirred cell fitted with 10kDa Millipore Ultrafiltration membrane.

### **SP- Sepharose Chromatography:**

15 Column matrix Q-Sepharose Fast Flow  
 Column dimensions 1.6 cm internal diameter x 15 cm bed height  
 Column volume 30 mL  
 Column preparation Packed in 0.1 M Sodium Hydroxide at 150 cm/hr  
 Operational flow rate 100 cm/hr (3.35mL/min)  
 20 Clean 0.1M Sodium Hydroxide (2 column volumes)  
 Loading capacity 10 mg /ml matrix  
 Equilibration 25mM Sodium acetate pH 5.00 /25mM NaCl (8 column volumes)  
 Post load wash 25mM Sodium acetate pH 5.00/25mM NaCl (6column volumes)  
 Elution 25mM Sodium acetate pH 5.00/186mM NaCl (25column volumes)  
 25 Strip buffer 25mM Sodium acetate pH 5.00/2M NaCl (2column volumes)  
 Wash 0.1M Sodium Hydroxide (2 column volumes)  
 30 24ml of buffer exchanged Q-Sepharose eluate was loaded onto the SP-Sepharose column at a loading capacity of 10mg/ml matrix. The column was operated in a binding mode; the eluate was collected as fractions. Fractions across the elution profile were analysed by GP-HPLC to determine aggregate levels results and are shown in Table 8. Samples following each chromatography step were collected and analysed for rProtein A residues, results are presented in Table 9.

**Table 8: rProtein A ELISA Results after each chromatography step**

Sample ID	rProtein A levels(ng/mg)	Antibody concentration (mg/ml)
MabSelect Protein A eluate (after conc/diaf)	2.64	46.7
Q-Sepharose eluate	<4	8.10
SP-Sepharose eluate pool F(1-16)	<4	1.79

5

**Table 9: GP-HPLC Analysis of SP-Sepharose fractions**

Sample ID	% Aggregates	Absorbance (A <sub>280</sub> )
SP eluate pool F(1-3)	0.57	11.2
SP eluate pool F(4-6)	1.10	2.7
SP eluate pool F(7-9)	2.07	0.655
SP eluate pool F(10-12)	2.12	0.351
SP eluate pool F(13-15)	2.56	0.208


10

Conclusion: Fractionation of aggregates was observed during step elution of the antibody peak; see Table 4. With the aggregate enriched fractions eluting later (on the tail end of the elution peak) compared to non-aggregate containing fractions. The tail fractions can be omitted from the main pool to obtain 99% monomeric pool and still have high recovery (>95%).

15

**Table 4 Fractionation of aggregates across the SP-Sepharose elution peak for high pI antibody**

20

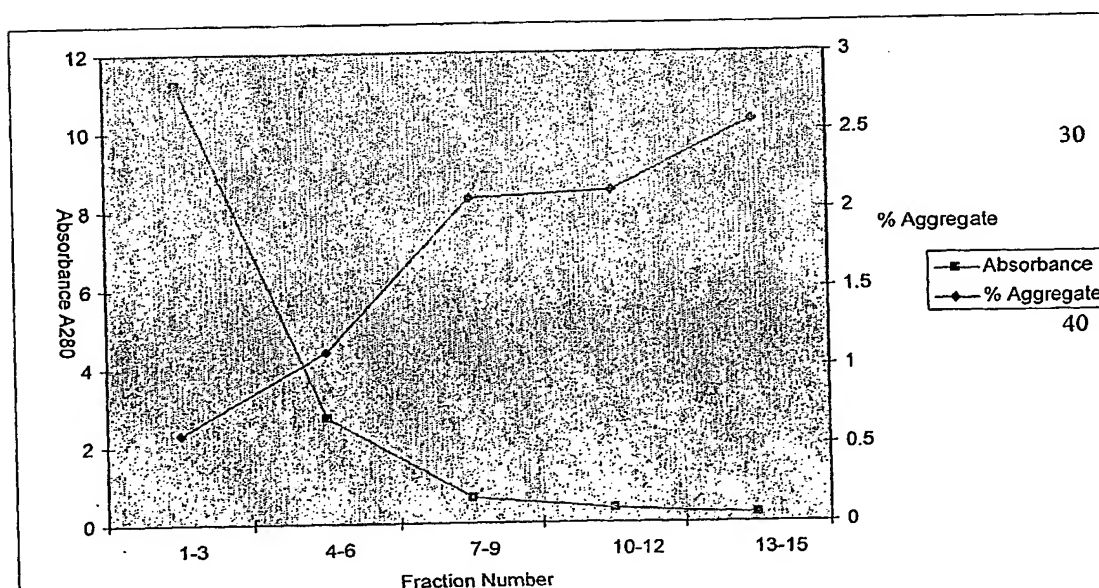

30

40

## Claims

1. Method of purifying an antibody, preferably an IgG antibody, comprising the steps of:
  - 5           • Purifying an antibody by means of protein A affinity chromatography wherein the protein A is a native protein A or a functional derivative thereof,
  - Loading the purified antibody on an ion exchange material under conditions which allow for binding of the protein A or its derivative and
  - 10          • Collecting the antibody, preferably collecting at least 70%, more preferably collecting at least 80%, most preferably collecting at least 90% of the amount of antibody loaded onto the ion exchange material in the flow-through of the ion exchanger whilst a contaminant protein A is bound to the ion exchange material
- 15           2. Method according to claim 2, characterized in that the protein A is a recombinant protein A that is engineered such as to allow of single-point attachment to a column material.
- 20           3. Method according to claim 2, characterized in that the recombinant protein A comprises a cysteine in its amino acid sequence.
4. Method according to claim 3, characterized in that the cysteine is comprised in a segment that consists of the last 30 Amino acids of the C-terminus of the amino acid
- 25           sequence of the recombinant protein A.
5. Method according to claim 3, characterized in that the recombinant protein A is attached by at least 50 % via a thioether sulphur bond to the chromatographic support material of the protein A affinity chromatography.
- 30           6. Method according to one of the preceding claims, characterized in that the protein A or its functional derivative is reduced to a concentration of <1ng/mg IgG in the flow-through of the ion-exchanger.

7. Method according to one of the preceding claims, characterized in that the ion exchanger is an anion exchanger.
8. Method according to one of the preceding claims, characterized in that the antibody has  
5 a pI of at least 6.5 or above.
9. Method of purifying an antibody comprising the steps of:
- Purifying an antibody by means of protein A affinity chromatography wherein the protein A is a native protein A or a functional derivative thereof,
  - 10 • Loading the purified antibody on a first ion exchange material under conditions which allow for binding of the protein A or its derivative,
  - Collecting the antibody, preferably collecting at least 70%, more preferably collecting at least 80%, most preferably collecting at least 90% of the amount of antibody loaded onto the ion exchange material in the flow-through of the ion  
15 exchanger whilst a contaminant protein A is bound to the ion exchange material,
  - And further purifying the antibody by loading on, binding to and eluting it from a second ion exchanger.
10. Method according to claim 9, ~~characterized in that the first ion exchanger is an anion~~  
20 ~~exchanger.~~
11. Method according to one claim 10, characterized in that the antibody has a pI of at least 7.5 or above.
- 25 12. Method according to claim 10, characterized in that the second ion exchanger is a cation exchanger.
13. Method according to claim 9, characterized in that the purified antibody is essentially monomeric, non-aggregated antibody.
- 30 14. Method according to one of the preceding claims, characterized in that the antibody is a monoclonal antibody, preferably an IgG antibody wherein the IgG antibody may be chimeric or CDR-grafted IgG antibody.

15. Method according to claim 14, characterized in that the IgG antibody is at least in those portions of the antibody that are relevant for binding to a protein A, preferably in its Fc portion that is relevant for binding to protein A, more preferably in the C $\gamma$ 2-C $\gamma$ 3 interface region of IgG comprising the binding sites for protein A, of such species origin and IgG subclass origin which origins allow for high affinity binding to protein A.
16. Method according to claim 15, characterized in that the IgG antibody is selected from the group comprising human IgG1, IgG2 and IgG4 with regard to the Fc portion of the antibody.
17. Method according to claim 14, characterized in that the antibody is harvested from a cell culture prior to purifying the antibody by means of protein A affinity chromatography.
18. Method according to claim 17, characterized in that the antibody is harvested from a mammalian cell culture.
19. Method according to claim 17 or 18, characterized in that the antibody that is to be purified by means of protein A affinity chromatography is not treated as to inactivate proteases, preferably is not in admixture with at least one protease inhibitor.
20. Method according to claim 19, characterized in that the protease inhibitor is selected from the group consisting of PMSF, a proteinase inhibiting peptide,  $\epsilon$ -caproic acid, and a reducing sulfhydryl compound.

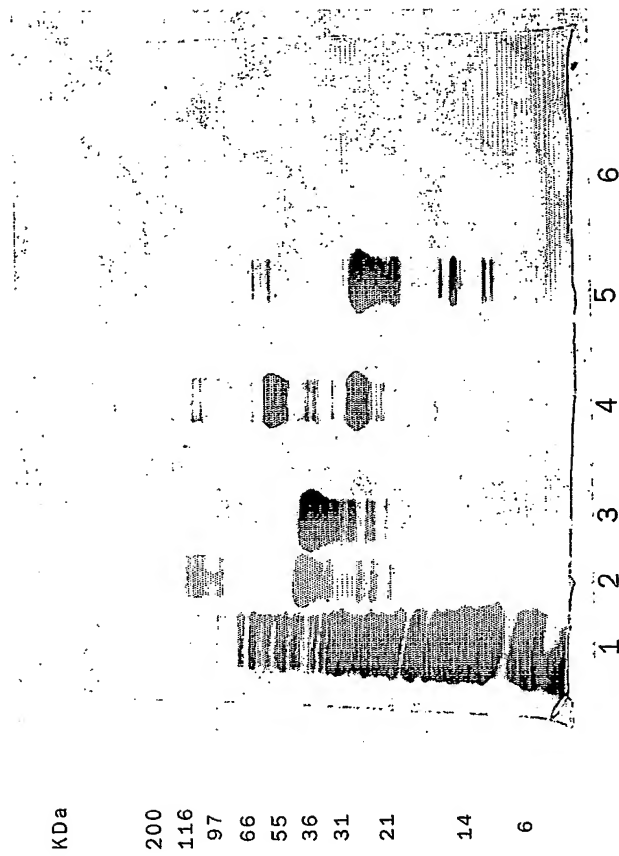
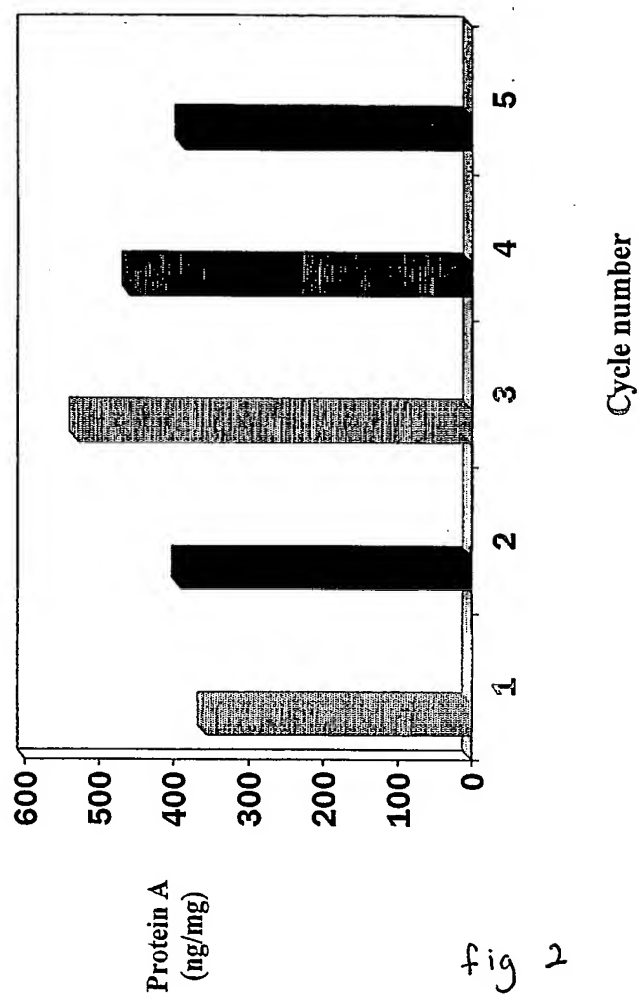


fig 1



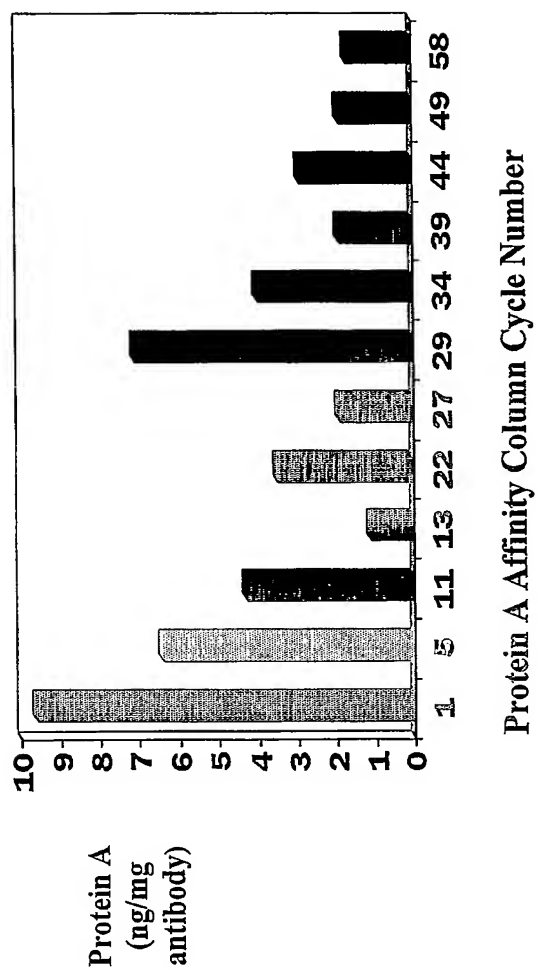


fig 3



# INTERNATIONAL SEARCH REPORT

International Application No  
PCT/EP2004/002041

A. CLASSIFICATION OF SUBJECT MATTER  
IPC 7 C07K16/06 C07K1/36  
//C07K17/02

According to International Patent Classification (IPC) or to both national classification and IPC

## B. FIELDS SEARCHED

Minimum documentation searched (classification system followed by classification symbols)  
IPC 7 C07K

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the International search (name of data base and, where practical, search terms used)

EPO-Internal, WPI Data, PAJ, BIOSIS

## C. DOCUMENTS CONSIDERED TO BE RELEVANT

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
X	US 4 983 722 A (BLOOM JAMES W ET AL) 8 January 1991 (1991-01-08) cited in the application column 1, line 54 - column 2, line 14 claims; figures; examples; tables	1,7,9, 10,13-15
X	US 5 429 746 A (ERICKSON JOHN C ET AL) 4 July 1995 (1995-07-04) column 3, line 37 - line 68 claims; examples	1,9, 13-20
X	JISKOOT W, ET AL.: "A purification strategy for the production of clinical grade monoclonal antibodies" DEVELOPMENTS IN BIOLOGICAL STANDARIZATION, vol. 71, 1990, pages 73-78, XP008031347 SWITZERLAND the whole document	1,7-11, 13-20
	----- -/-	

☒ Further documents are listed in the continuation of box C.

☒ Patent family members are listed in annex.

### \* Special categories of cited documents:

- \*A\* document defining the general state of the art which is not considered to be of particular relevance
- \*E\* earlier document but published on or after the international filing date
- \*L\* document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)
- \*O\* document referring to an oral disclosure, use, exhibition or other means
- \*P\* document published prior to the international filing date but later than the priority date claimed

- \*T\* later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention
- \*X\* document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone
- \*Y\* document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art.
- \*G\* document member of the same patent family

Date of the actual completion of the international search

22 June 2004

Date of mailing of the international search report

05/07/2004

Name and mailing address of the ISA

European Patent Office, P.B. 5818 Patentlaan 2  
NL - 2280 HV Rijswijk  
Tel. (+31-70) 340-2040, Tx. 31 651 epo nl,  
Fax: (+31-70) 340-3016

Authorized officer

Rankin, R

# INTERNATIONAL SEARCH REPORT

International Application No  
PCT/EP2004/002041

C.(Continuation) DOCUMENTS CONSIDERED TO BE RELEVANT		
Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
X	BLANK GS, ET AL.: "Expanded bed absorption in the purification of monoclonal antibodies: a comparison of process alternatives." BIOSEPARATION, vol. 10, no. 1-3, 2001, pages 65-71, XP008031348 NETHERLANDS the whole document	1,6-9, 13-20
A	----- DATABASE WPI Section Ch, Week 199533 Derwent Publications Ltd., London, GB; Class B04, AN 1995-250742 XP002283061 & JP 07 155194 A (TEIJIN LTD) 20 June 1995 (1995-06-20) abstract	1-20
A	----- EP 0 284 368 A (REPLIGEN CORP) 28 September 1988 (1988-09-28) cited in the application page 2, line 55 - page 3, line 5 examples 8-10 claims 10-12	1-20
A	----- US 6 399 750 B1 (JOHANSSON INGEMAR) 4 June 2002 (2002-06-04) column 1, line 15 - line 21 examples	1-20
	-----	

# INTERNATIONAL SEARCH REPORT

Information on patent family members

International Application No

PCT/EP2004/002041

Patent document cited in search report		Publication date	Patent family member(s)	Publication date
US 4983722	A	08-01-1991	AT 121752 T CA 1339188 C DE 68922338 D1 DE 68922338 T2 EP 0345549 A2 ES 2071628 T3 JP 2084193 A JP 2664478 B2 US 5115101 A	15-05-1995 29-07-1997 01-06-1995 31-08-1995 13-12-1989 01-07-1995 26-03-1990 15-10-1997 19-05-1992
US 5429746	A	04-07-1995	AT 218387 T AU 689552 B2 AU 1843395 A BR 9507100 A CA 2183888 A1 CN 1146730 A CZ 9602481 A3 DE 69526929 D1 DE 69526929 T2 DK 746398 T3 EP 0746398 A1 ES 2177632 T3 HK 1013806 A1 HU 74845 A2 JP 9509658 T NO 963475 A NZ 281480 A PT 746398 T SI 746398 T1 WO 9522389 A1 ZA 9501372 A	15-06-2002 02-04-1998 04-09-1995 16-09-1997 24-08-1995 02-04-1997 16-04-1997 11-07-2002 02-01-2003 16-09-2002 11-12-1996 16-12-2002 07-02-2003 28-02-1997 30-09-1997 21-10-1996 26-06-1998 31-10-2002 31-10-2002 24-08-1995 24-10-1995
JP 7155194	A	20-06-1995	NONE	
EP 0284368	A	28-09-1988	US 5084559 A EP 0284368 A1 JP 63267281 A	28-01-1992 28-09-1988 04-11-1988
US 6399750	B1	04-06-2002	AU 709290 B2 AU 7593096 A CA 2236900 A1 DE 69629127 D1 DE 69629127 T2 EP 0873353 A1 JP 2000500649 T WO 9717361 A1	26-08-1999 29-05-1997 15-05-1997 21-08-2003 15-04-2004 28-10-1998 25-01-2000 15-05-1997

**THIS PAGE BLANK (USPTO)**

**This Page is Inserted by IFW Indexing and Scanning  
Operations and is not part of the Official Record**

**BEST AVAILABLE IMAGES**

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☐ BLACK BORDERS
- ☐ IMAGE CUT OFF AT TOP, BOTTOM OR SIDES
- ☐ FADED TEXT OR DRAWING
- ☐ BLURRED OR ILLEGIBLE TEXT OR DRAWING
- ☐ SKEWED/SLANTED IMAGES
- ☐ COLOR OR BLACK AND WHITE PHOTOGRAPHS
- ☐ GRAY SCALE DOCUMENTS
- ☒ LINES OR MARKS ON ORIGINAL DOCUMENT
- ☐ REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY
- ☐ OTHER: \_\_\_\_\_

**IMAGES ARE BEST AVAILABLE COPY.**

**As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.**

**THIS PAGE BLANK (USPTO)**